

Vacuum Assisted Delivery System




RELIABILITY AND FLEXIBILITY



Medela Basic for Vacuum Assisted Delivery

Reliable vacuum assisted delivery

The Medela Vacuum Delivery System offers reliability and flexibility whenever the birth process requires vacuum assistance. Together with the Medela Basic suction pump and your choice of rigid or soft, reusable or disposable suction cups and a wide range of accessories, the system delivers effective performance. The Basic suction pump offers fast and efficient vacuum build-up, easy maintenance and precise vacuum settings. By producing a consistent active vacuum to keep the suction cup in place, it also may reduce the risk of cup pop-offs.

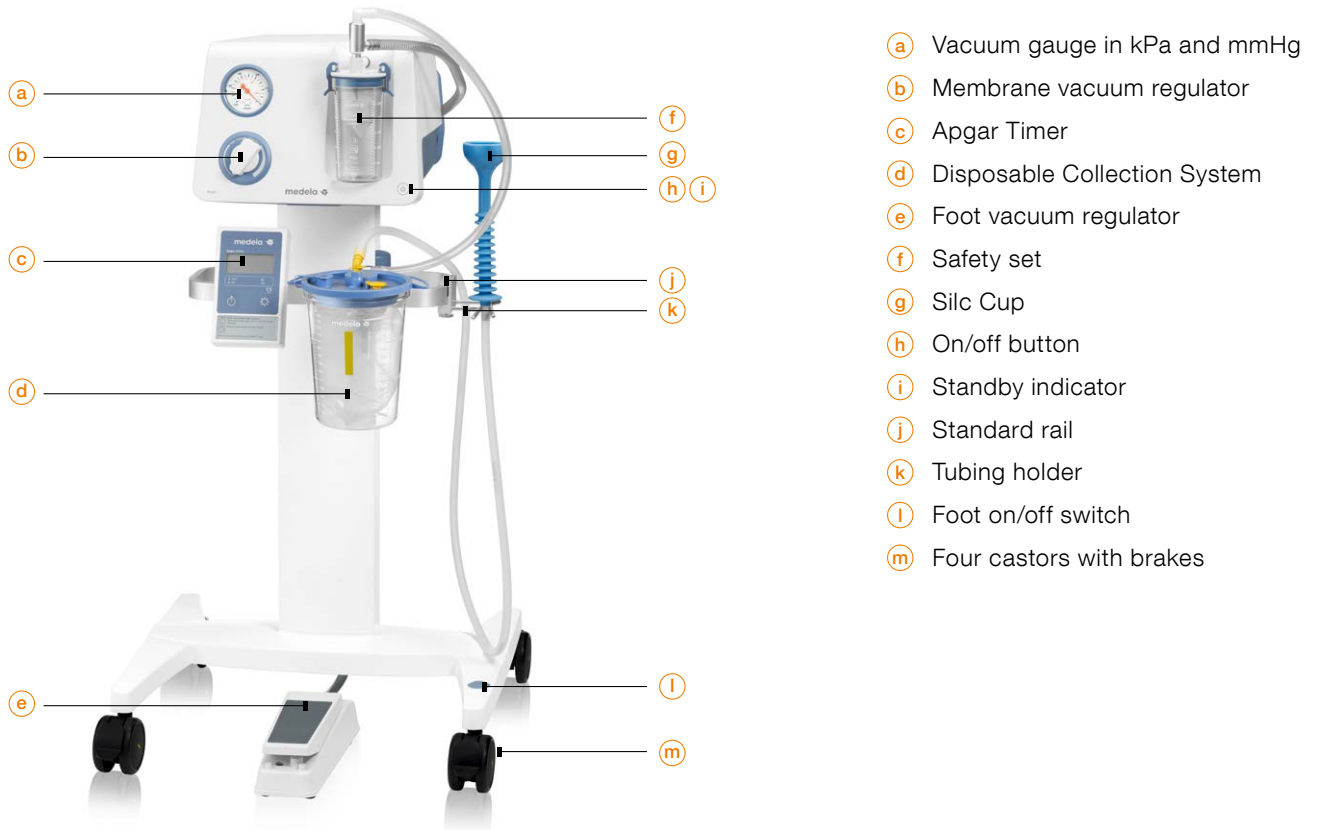

- a Vacuum gauge in kPa and mmHg
- b Membrane vacuum regulator
- c Apgar Timer
- d Disposable Collection System
- e Foot vacuum regulator
- f Safety set
- g Silc Cup
- h On/off button
- i Standby indicator
- j Standard rail
- k Tubing holder
- l Foot on/off switch
- m Four castors with brakes

Fast Facts

- Medela offers healthcare professionals 50 years of expertise in medical vacuum technology.
- Our swiss development and manufacturing process yields mechanical precision to guarantee excellent performance.
- Electric vacuum generation at the touch of a button ensures fast vacuum build-up.
- Active vacuum ensures firm attachment on the fetal head, potentially avoiding inadvertent pop-offs.
- Quick setup since the compatible components fit together securely.
- Choice of type of suction cups and size to ensure gentle and individual vacuum assisted delivery.


Medela's suction cup product range:

- Extensive range of suction cups to fit specific needs
- Cups and tubing compatible with Basic suction pump and other manual or electric suction pumps

Rigid Suction Cups – for optimal adhesion

Malmström Cups

- Proven effective for decades
- Used when the fetus is in the occipito-anterior position
- Cup design allows for good seal to help minimise pop-offs
- Autoclavable stainless steel
- Available in 40, 50 and 60 mm sizes

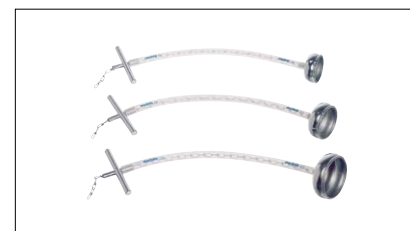

Malmström Cups

Bird Cups

- Most widely sold suction cup
- Based on the design of the Malmström Cup, with separate suction and traction ports
- Autoclavable stainless steel
- Available in 40, 50 and 60 mm sizes



Bird Cups

Bird Posterior Cup


- Especially designed for use when the fetus is in the occipito-posterior position
- Vacuum port on side of suction cup and low profile to ensure easier insertion
- Autoclavable stainless steel
- 50 mm size



Bird Posterior Cup

Single Use Bird Cup

- Polycarbonate version of the reusable Bird Cup
- Traction handle can be tilted for use when the fetus is in the occipito-anterior or posterior position
- 50 mm size

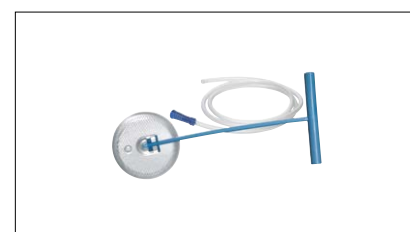


Single Use Bird Cup

Soft Suction Cups – for gentle vacuum assisted delivery

Reusable and Single Use Silc Cups

- Attach gently to fetal scalp to minimize risk of fetal head trauma
- One-piece design for easy handling and use
- Smooth surface with longitudinal ridge makes it possible to observe rotation
- Reusable cups made of medical grade silicone, available in 50 and 60 mm sizes
- Single Use Cup made of thermoplastic elastomer, size 60 mm

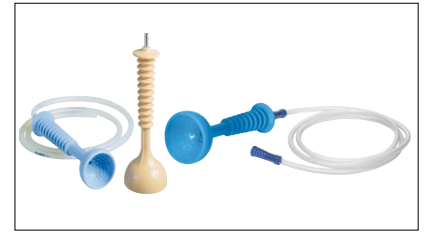

Silc Cups

CaesAid Cup for Caesarean sections

- Gentle support for delivering fetal head during Caesarean sections
- Can be attached to the pump used for suctioning during the operation
- Autoclavable medical grade silicone
- 60 mm size



CaesAid Cup

Special accessories

- Foot vacuum regulator for hands-free increase and quick release of vacuum
- Reusable or Disposable Collection System
- Apgar Timer
- Holders for tubing and accessories

Medela offers healthcare professionals 50 years of experience in medical vacuum technology. Our Swiss development and manufacturing process yields mechanical precision to guarantee excellent performance.

Technical data – Basic suction pump



high vacuum
– 90 kPa
– 675 mmHg
Tolerance: +/- 15 %



H x W x D (rack version)
210 x 305 x 375 mm



ISO 9001
ISO 13485
CE (93/42/EEC), IIa

CE 0123



high flow
30 l/min.
+/- 5 l/min.



Power
100 – 240V



5 years warranty



9.3 kg (rack version)



IP21

CAUTION: Federal US law restricts this device to sale by or on the order of a physician.

Medical Vacuum Technology for Healthcare Professionals

Please contact us or your local Medela representative for details.

 Medela AG
Lättichstrasse 4b
6341 Baar, Switzerland
www.medela.com

Local contact:

Australia
Medela Pty Ltd
Medical Technology
3 Arco Lane
Heatherton, Vic 3202
Australia
Phone +61 3 9552 8600
Fax +61 3 9552 8699
contact@medela.com.au
www.medela.com.au

Canada
Medela Inc.
4160 Sladeview Cres., Unit #8
Mississauga, Ontario
L5L 0A1
Canada
Toll Free 800 435 8316
Toll Free Fax 800 995 7867
info@medela.ca
www.medela.ca

India
Medela India private limited
c/o Vatika Business Park
First floor, tower 2, Sohna
Road, Sec-49,
Gurgaon 122 002
India
Direct +91 124 4416999
Fax +91 124 4416990
www.medela.in

UK
Medela UK Ltd.
Huntsman Drive
Northbank Industrial Park
Irlam, Manchester M44 5EG
UK
Phone +44 161 776 0400
Fax +44 161 776 0444
info@medela.co.uk
www.medela.co.uk

USA
Medela Inc.
1101 Corporate Drive
McHenry, IL 60050
USA
Phone +1 877 735 1626
Fax +1 815 307 8942
info-healthcare@medela.com
www.medela-healthcare.us